



MEDTECH®
Leader in Health & Home care!

OXY GARD OG-01



PULSE OXYMETER


Features:

- Light & Compact design
- Quick measure
- Easy Operate & Carry
- Automatically power off
- Low Power Consumption
- Dual color OLED display


Useful for:

- Monitor Oxygen Saturation of Patient & during sports & high altitude activities.
- Monitor Heart rate
- Managing Asthama & COPD
- Managing supplemental Oxygen Use.

How to Use:

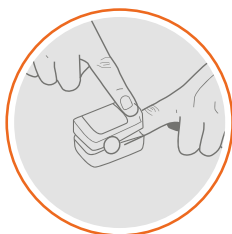

Place one of your fingers into the pulse oxymeter to the end and press the switch button one time on front panel to turn it on.

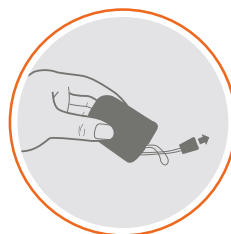

Tie the lanyard through the hole in the rear of the pulse oxymeter



OXY GARD OG-01 PULSE OXYMETER

Specifications:

- **Display Type** : OLED display
- **SpO2**
 - Display range : 0%~100%
 - Measurement range : 70%~100%
 - Accuracy : 70%~100%±2%; 0%~69% no definition
 - Resolution : 1%
- **Pulse Rate**
 - Display range : 30bpm 250bpm
 - Measure range : 30bpm 250bpm
 - Accuracy : 30bpm 99bpm, ±2bpm; 100bpm 250bpm, ±2%
 - Resolution : 1bpm
- **Perfusion Index**
 - Display range : 0.1% 20.0%
 - Measure range : 0.2% 20.0%
 - Accuracy : 0.2% 1.0%, ±0.2digits; 1.1% 20.0%, ±20%
 - Resolution : 0.1%

Probe Led Specifications:


	WAVELENGTH	RADIANT POWER
RED	660±3nm	3.2mW
IR	905±10nm	2.4mW



Pulse Rate



Oxygen Saturation

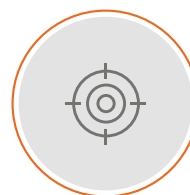

High Accuracy

Note: The information provided is best to our knowledge and subject to change without notice
For more Information, please contact:

MEDTECH LIFE PVT. LTD.

(Formerly Nulife Global Medical Devices Pvt. Ltd.)

B6, Byculla Service Industries,

D.K. Marg, Byculla, Mumbai 400 027, India.

T: +91 22 23748372/73 / +9122 6657 8989

E: info@medtechlife.com W : www.medtechlife.com



An ISO 9001: 2008
Certified Company