



Hand Rehabilitation Robot

SY-HRC12



Soft robotic
exoskeleton technology



6 training modes, cover
all stages of hand rehabilitation




Stronger power,
9 levels intensity to adjust



Core advantages

Combined with soft robotic exoskeleton technology and neuroscience, syrebo[®] hand rehabilitation robot is fully upgraded from function, power, gloves, display screen and so on, so as to improve hand mobility and accelerate the process of hand function rehabilitation from three levels of nerves, brain and muscles.

Suitable for patients with hand dysfunction caused by stroke, brain injury, hand trauma, cerebral palsy, orthopedic surgery, lymphedema after breast cancer treatment , etc.

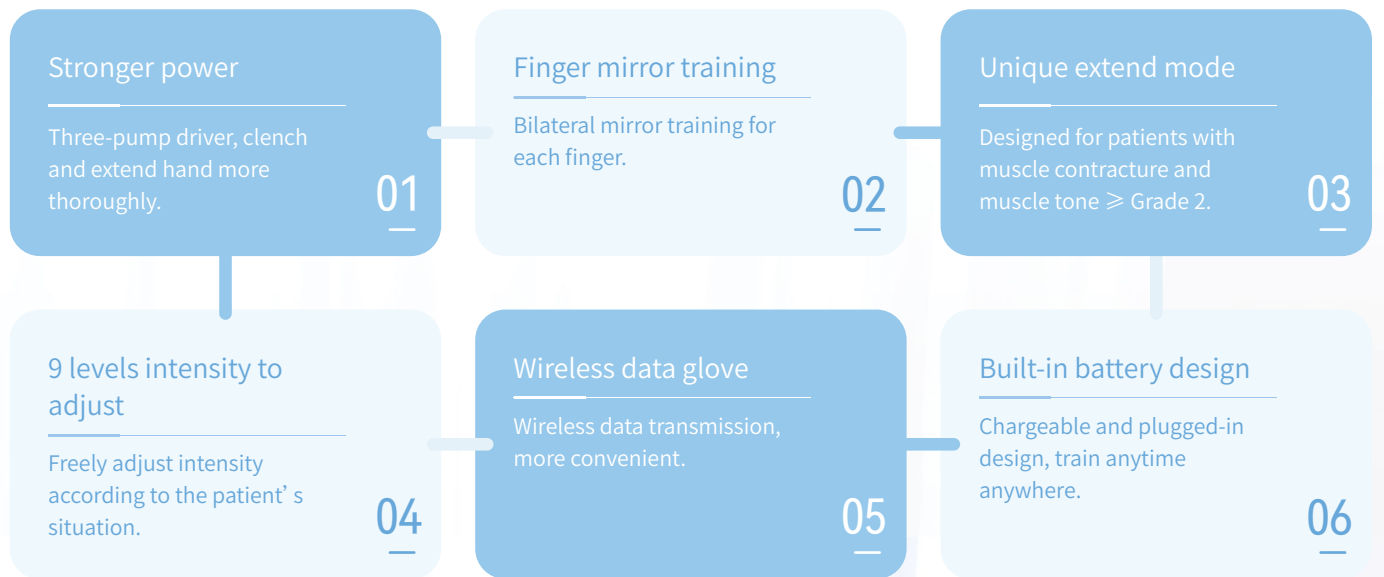

6 training modes



Passive Training

Power glove drives the affected hand to execute flexion and extension, and the glove can intelligently adjust the speed of flexion and extension.



Bilateral Mirror training

The healthy hand leads the affected side fingers to do synchronized actions.



Pinch Training

The thumb opposes the index finger, middle finger, ring finger and little finger in turn.



Extend Training

Static stretching hand to decrease muscle tone and increase range of motion.



Fine Motor Training

Power glove drives each affected finger to perform the flexion and extension.



Functional Training

Two or more fingers pinch training to simulate daily item pinch.

Product specifications

| Product name | Model | Size | Rated charging voltage | Rated charging current |
|-----------------------------------|----------|----------------|------------------------|------------------------|
| Syrebo® Hand Rehabilitation Robot | SY-HRC12 | 135*145*65(mm) | 9V | 1.5A |

Packing list

| Name | User manual (including warranty card) | Host | Power glove | Data glove | Quality certificate | Adaptor |
|----------|--|------|-------------|------------|---------------------|---------|
| Quantity | 1 | 1 | 1 | 1 | 1 | 1 |

80⁺ Countries & regions **4000⁺** Medical units **60000⁺** Families

Providing innovative rehabilitation products & services
for hospitals/communities/families



LinkedIn



Facebook

Shanghai Siyi Intelligent Technology Co., Ltd.

Website: www.syrebo.com

E-mail: info@syrebo.com

Mobile&WhatsApp: (+86) 17317983724

Office: (+86) 400-6886-275

Address: 4F, Building 13, 2388 Xiupu Road, Pudong New Area, Shanghai